

HEALTHY, GLOWING, GLAMOROUS SKIN SHOULD NOT BE A LUXURY

Stress, pollution, UV radiation and an unbalanced diet can all have a negative impact on our skin. As we age, our skin becomes dull, wrinkles begin to appear and the texture of our skin can become coarse and uneven. Proper skin care is essential for maintaining the vitality and glow of healthy, youthful skin.

Alma Beauty, a sub-brand of medical technology leader Alma Lasers, has developed an advanced skin rejuvenation solution that treats skin imperfections with instant results. ALMA BEAUTY REBOOST.

ALMA BEAUTY REBOOST uses the power of acoustic oscillations to enhance the effectiveness of the entire skin treatment session- from cleansing and exfoliation to the application and penetration of therapeutic dermatological ingredients.

REBOOST emits acoustic waves of mechanical pressure on the skin. These create a **push and pull effect**, exfoliating the skin by removing dead skin cells and excess sebum. The acoustic vibrations also cause gaseous bubbles within the skin's barrier, or the stratum corneum,

to rapidly expand and collapse. This disorganizes the lipid bi-layer, creating crevices and pores between 25-75 micrometers in diameter, and increases the fluidity of the lipids in the stratum corneum.

More fluidity in sub-epidermal tissue allows topical ingredients to penetrate the skin more deeply and effectively. Both hydrophobic and hydrophilic compounds are effectively delivered to the target tissue depth and remain active for up to 30 minutes after application, for maximally effective skin repair and aesthetic enhancement.

Multiple clinical studies and extensive laboratory testing has proven the ability of REBOOST technology to deliver active ingredients deep beyond the skin barrier and has shown positive change in skin appearance following each treatment.

The REBOOST treatment is easy to administer, is safe for all skin types and can be performed all year round. Patients experience minimal to no discomfort or downtime and enjoy immediate, dramatic results.





Facial rejuvenation treatments are an absolute must. They are part of my routine skin care regimen" Alberta Ferretti, Fashion Designer


Curtsey of Novamedical, Clinical Department

REBOOST is now offered with the **PixelPen** and **Youth B(ox**) treatment kit, by Mimi Luzon, a skin care expert with 25 years of experience in the field..

PixelPen is a fractional vibrating micro-needling device that perforates the skin before treatment. These perforations, in combination with REBOOST's acoustic pressure technology, allows for maximum trans-epidermal delivery of cosmeceuticals and active ingredients to the deeper layers of the skin.

THE PIXELPEN FEATURES:

Universal, sterilized and disposable tip

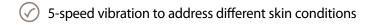

~2mm multiple stainless steel needles that can be pre-configured to a specific depth of penetration

The PixelPen is ideal for enhanced treatment of cutaneous indications such as skin tightening and rejuvenation, acne scars, fine lines and wrinkles and reducing pore size.

Delivering the benefits of the **Youth B(ox)** series has never been easier - with REBOOST.

ALMA BEAUTY REBOOST is now offered with the Youth B(ox) treatment kit by Mimi Luzon to optimally maintain youthful skin at any age.

Youth B(ox) features special ingredients* at the forefront of anti-aging skin treatments.

All of the formulas have been time-tested, yielding superior improvement of the skin's texture and tone within a relatively short period of time.

100% Pure Hyaluronic Acid

Cyber C Vitamin C Serum

HGF Youth Serum

Collagen Synthesizer



With Mimi, I feel that my skin care is in the best hands and that feeling cannot be replaced. Every treatment gives me a refreshed, rejuvenated complexion."

Eva Cavalli, Cavalli Business Partner



^{*} None of the products in the series includes parabens. None were tested on animals. All products have been dermatologically tested.



Delivering the benefits of the Youth B(ox) series has never been easier - with REBOOST:

Safe For Use All Year Round

Suitable For All Skin Types

Improve Skin Imperfections

Immediate Results

Minimal to No Downtime

ALMA BEAUTY REBOOST is ideal for a variety of indications including:

Skin Rejuvenation

Skin Whitening

Skin Tightening

Skin Imperfections:

• Oily or Problematic Skin

Skin • Dehydrated Skin

• Fine Lines and Wrinkles

Redness

• Sun Damage



Natural Beauty powered by Medical Technology

Alma Beauty combines the sensitivity and delicate approach of aesthetic care with revolutionary advancements from the field of aesthetic medicine, to offer breakthrough treatment solutions that achieve natural-looking beauty with aweinspiring results. Alma Beauty solutions have a superb safety record proven in extensive clinical studies, with safe, highly effective treatments for a variety of beauty needs.

Alma Beauty is a sub-brand of Alma Lasers, a global innovator of laser, light-based, radiofrequency and ultrasound technologies for the aesthetic and surgical markets.

Reveal the full potential of your natural beauty.
To learn more about Alma Beauty, visit us at:
www.almabeauty.com

©2015 Alma Lasers Inc. All rights reserved. Alma, Alma Lasers, the Alma Lasers logo, Alma Beauty, the Alma Beauty logo, REBOOST registered trademarks or trademarks of Alma Lasers Inc. in the United States and/or other countries.

www.almalasers.com

Alma Lasers AT GmbH

AUSTRIA, Peuerbachstraße 2, Top 208, 4040 Linz Phone: +43 (0) 732 / 70 05 07 Fax: +43 (0) 732 / 70 00 13 marketing@almabeauty.com info@alma-lasers.de www.alma-lasers.de



